

Labo Disk Holder

LDH Type

Major Applications

Lab-scale evaluation of various liquids



Features

- Compatible with 25mm diameter disc filter media (maximum thickness 5mm)
- Luer lock joint


Advantages

- Filters can be selected according to the purpose.
- Ideal for filtration of very small volumes
- Easy to attach / detach to / from syringe etc.

Picture of parts




Internal structure

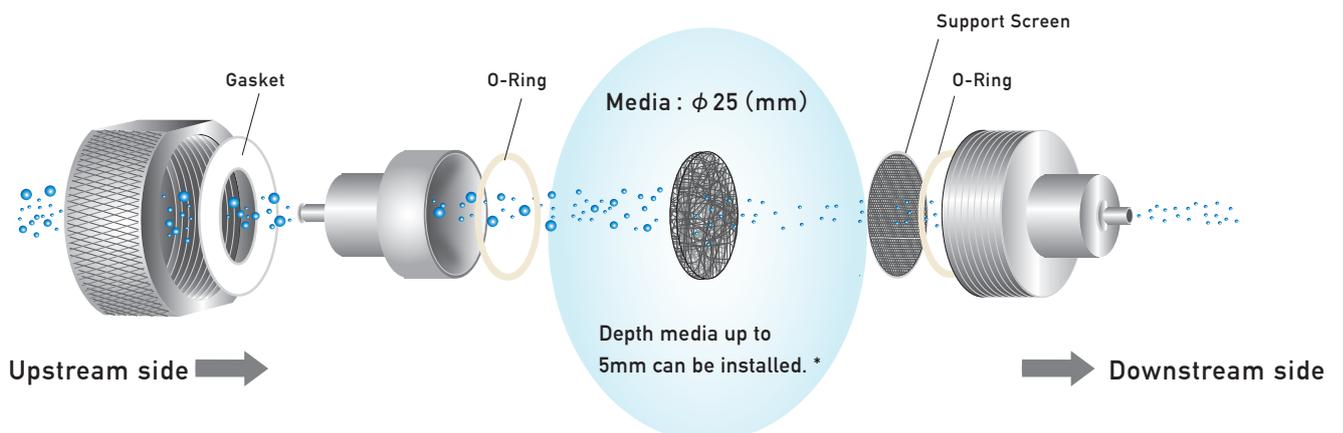

*Please refer to the disk filter catalog for the settable filter media.

Specifications

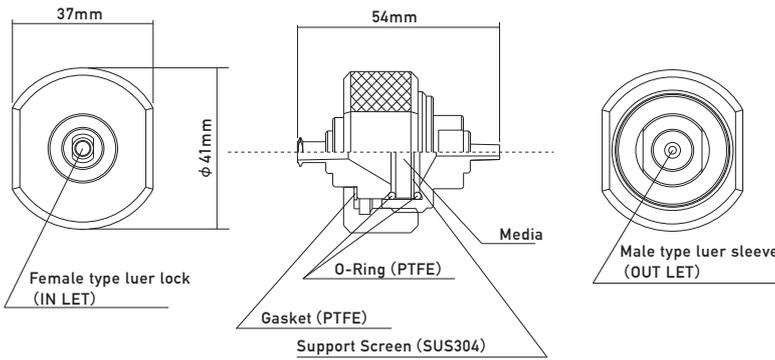

Disk Media	O.D. ϕ 25mm Media up to 5mm
Holder Material	SUS304
Connection diameter (IN LET)	Female type luer lock
Connection diameter (OUT LET)	Male type luer sleeve
Gasket	PTFE
O-Ring	PTFE
Support Screen	SUS304
Internal Volume of the Holder	7ml
Surface Treatment	Inside : Buff (#400)
Design Pressure	0.49MPa (71psi)
Design Temperature	40°C (104°F)

Ordering Information

LDH-L1

*The contents of the catalog are subject to change without notice.

*The performance data listed in the catalog are Typical values obtained under specific conditions based on our tests.

ROKI GROUP International Pte. Ltd.

6-20-12, Minami-Oi, Shinagawa-ku Tokyo, 140-0013 Japan
TEL: +81-3-5764-1131 FAX: +81-3-5764-0681

www.rokiglobal.com

Manufacturing is based on our Quality Management Systems that meet ISO9001 standards.

Scope

Design, Development, manufacture, and sales of filter cartridges, housings and filtration equipment.



JQA-QMA16323

MS CM009

16th Issue
LDH200702E